

# Druglizer

LE5 drug tester



## THE ULTIMATE HANDHELD DRUG TESTING SOLUTION

- ✓ Rapid, accurate results from 90 seconds
- ✓ Electronic cartridge batch validation
- ✓ Clear, visual electronic display
- ✓ Cost-effective drug testing
- ✓ Easy to use operation

Omega Laboratories is proud to be the Master Distributor for Alcolizer Technology products. The Druglizer LE5 allows you to add a rapid oral fluid drug detection device to your current hair testing program.

### ACCURATE, FAST, LAW ENFORCEMENT

#### GRADE DRUG TESTING

LE5 Drug Tester technology incorporates state-of-the-art Linear Detector Arrays and a patented Optical Filtering System to deliver accurate and fast drug test results from just 90 seconds. Results are clearly defined and displayed so the user is not required to interpret the presence of a line like other test devices.

#### MULTI-DRUG GROUPS TESTING

Developed to detect a variety of drug groups, LE5 Drug Tester utilizes lateral flow to provide accurate instant POCT results.

#### FAST, EASY CALIBRATION

Maximum unit availability – the module exchange system enables easy on-site calibration. Achieve long term savings via reduced unit management and no operational disruption.

### SERVICE AND SUPPORT

Each LE5 Drug Tester module has a product lifetime warranty on the parts when maintained under an AlcoCARE™ plan. The LE5 Drug Tester is supported by a global network of service providers and factory trained professional technicians for service and calibration every 6 months.

### ROBUST AND DURABLE

Designed and made for demanding conditions, the LE5 Drug Tester delivers outstanding performance in rugged and remote environments.



# Druglizer

LE5 drug tester



## DRUGLIZER LE5 TECHNICAL SPECIFICATIONS

**INSTRUMENT APPLICATION** ▶ Law enforcement and Industry

**SENSOR TYPE** ▶ Lateral flow strips with patented Optoelectronic sensor

**ACCURACY** ▶ Designed to Australian Standards AS 4760-2006 (Procedures for specimen collection and the detection and quantitation of drugs in oral fluid)

**TEST TIME** ▶ Test results from 90 seconds

**DISPLAY** ▶ Transreflective colour graphics LCD screen

**MOUTHPIECE** ▶ Druglizer Oral Fluid Collector

**CALIBRATION** ▶ Module exchange system for on-site calibration

**POWER SUPPLY** ▶ Rechargeable Li-Ion battery system

**WEIGHT** ▶ 310g (including battery and cover)

**DIMENSIONS** ▶ 196mm H x 66mm W x 43mm D

## LE5 Easy On-site Sample Module Exchange Calibration



Every 6 months, Alcolizer will send a calibrated Druglizer LE5 replacement




Using the supplied allen key, remove the screw from the back of the LE5.



Slide off the expired module.



Slide on the calibrated module. Replace the screw. Turn your Druglizer LE5 on – it is now ready to use!

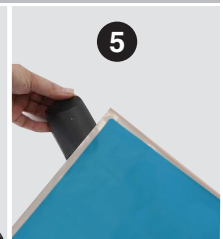

Return the expired module back in the supplied post bag.